



LorX[®]

cervical peek disc prosthesis

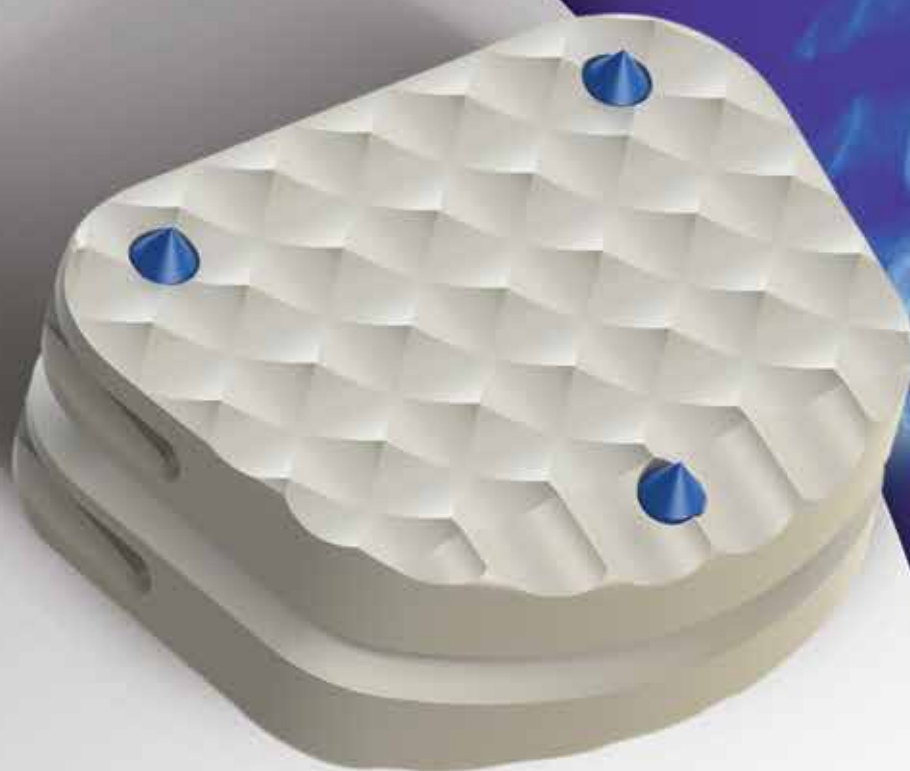

 Tria Spine[®]
SPINE | ORTHOPEDICS


DESIGN FEATURES

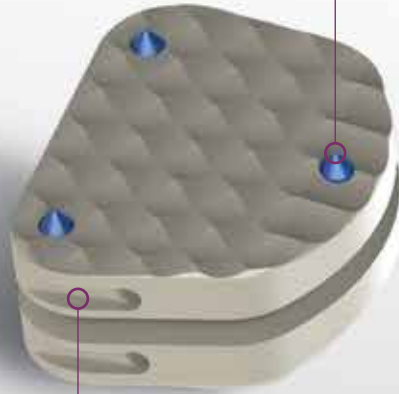

Titanium pins
specially designed six titanium pins for perfect fixation


Better Stability
two instrument sockets for extra stability during insertion

Anatomical Shape
original and anatomical design for a secure and better placement

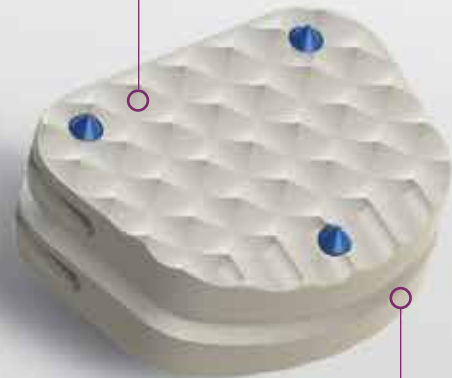

Movement
special angulation mechanism which allows 10° both ways in flexion - extension, 7.5° both ways in lateral bending

Biocompatibility
radiolucent and biocompatible PEEK material (VESTAKEEP by Evonik, ASTM F2026) which provides modulus of elasticity similar to the bone

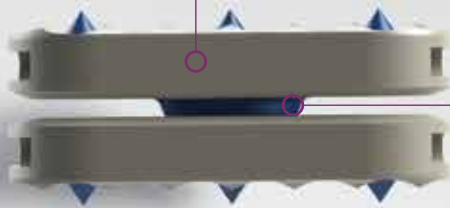

UHMWPE core
enhanced core which allows flexion, extension and rotation

LorX[®] Cervical PEEK Disc Prosthesis

LX-10-525	15 x 12 mm	5.5 mm
LX-10-526	15 x 12 mm	6.5 mm
LX-10-527	15 x 12 mm	7.5 mm